



Programming of the 7955 times used in Oxford Cryo-Pump Controllers

Note: Programming to be done with No Power to Rack. For **Timer 1** the desired configuration is:

No Power to Rack

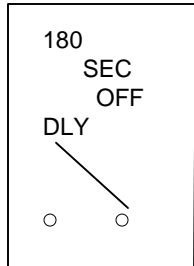


With Rack Active

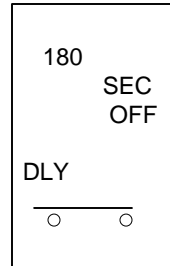

Goes to this
configuration
when Rack
Active


A. To start programming mode

Press “prog” button and hold for 3-5 seconds. At the end of that time display will show.

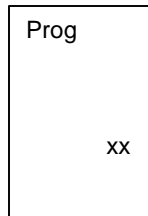

B. To set timer relay configuration


Press “Prog” button as many times as necessary to display.

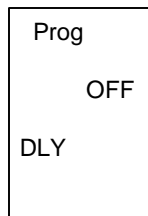

Press “Reset” button when you have the desired relay configuration.

C. To set time scale (hrs, min, sec):


Press "Prog" button as many times as necessary to display.

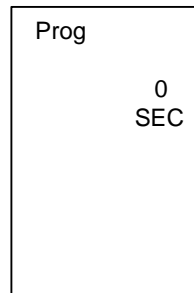

Press reset button when you have the desired time scale


D. To set delay time (least significant digit)

- The desired setting is 180, least significant digit can be set by pressing "Prog" button until desired # is set, when set press "Reset" button and display will move to next significant digit

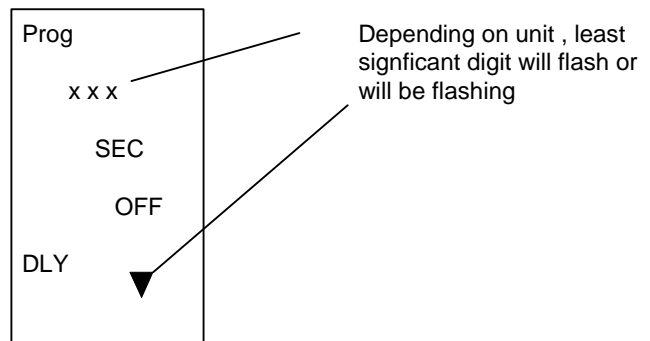


E. To set delay time (next significant digit)

- Press "Prog" button until the desired # is set, when set press "Reset" button and the display will move to the next significant digit

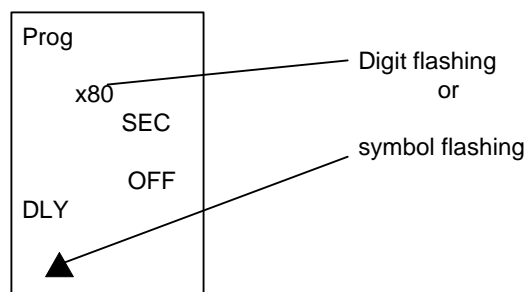


F. To set delay time (most significant digit)

- Press "Prog" button until the desired # is set
- Press "Reset" button

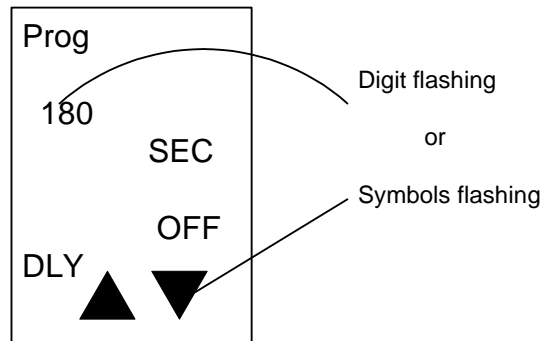

G. Set trigger type [Level (LEU) or Edge (Edg)]


- Press "Prog" button to display LEU
- Press "Reset" button when desired trigger type is indicated

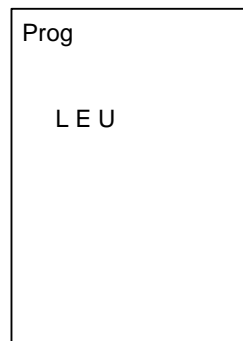

H. The timer now reviews your selected setting (if changes were made

{Timer Relay Conf.}


If the selection is correct
- Press "Reset" button

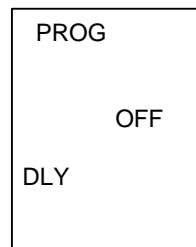

If incorrect press "Prog" button until you have the correct selection, then press "Reset"

{Time Scale}


If selection is correct
- Press "Reset" button

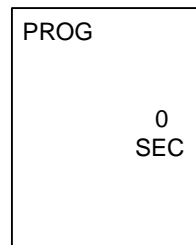

If incorrect press "Prog" button until you have the correct selection then press "Reset" {time setting (Isd)}


{Time setting Isd}

If selection is correct
- Press "Reset" button

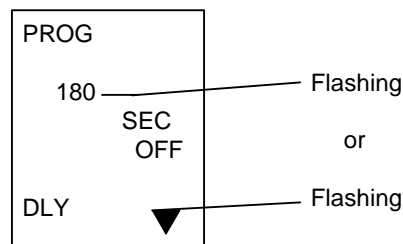

If incorrect press "Prog" button until you have the correct selection, then press "Reset"


If selection is correct
- Press "Reset" button

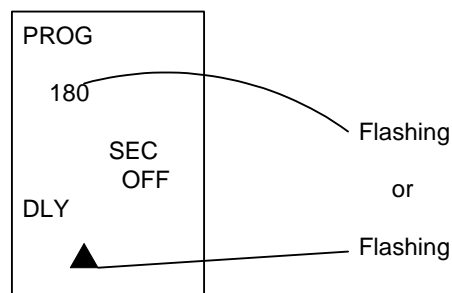

{Time setting nsd}

If incorrect press "Prog" button until you have the correct selection - then press "Reset"


If selection is correct
- Press "Reset" button

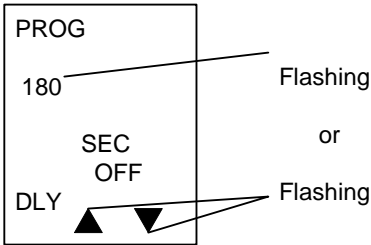

[Time Setting MSD]

If incorrect, press "Prog" button until you have the correct selection, then press "Reset"

If selection is correct
- Press "Reset" button

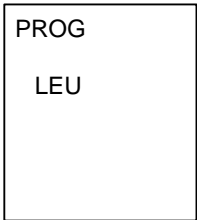


[Trigger type]


If incorrect, press "Prog" button until you have the correct selection, then press "Reset"

Operational state


No Rack Active

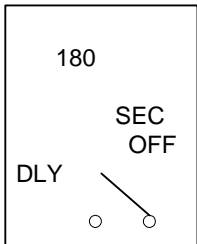


Timer 2 - The same sequence is now done for Timer 2. The desired configuration is [Note - Trigger type is Edg], and desired configuration is the same Rack Active or not.

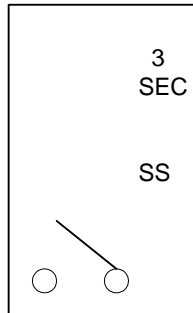

- A. To start programming mode:
Press "Prog" button and hold for 3-5 seconds, at the end of that time display will show

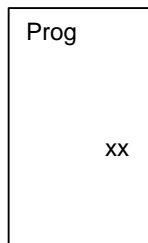


- B. To set timer relay configuration:
Press "Prog" button as many times as necessary to display

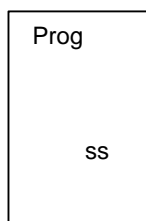

Press "Reset" button when you have the desired relay configuration.


- C. To set time scale (hrs, min, sec):
Press "prog" button as many times necessary to display

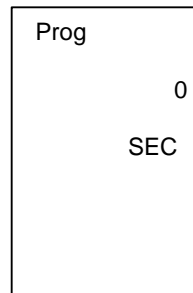

Press reset button when you have the desired time scale.


- D. To set delay time (least significant digit)

- The desired setting is 003, least significant digit can be set by pressing "Prog" button until desired # is set, when set press "Reset" button and display will move to next significant digit

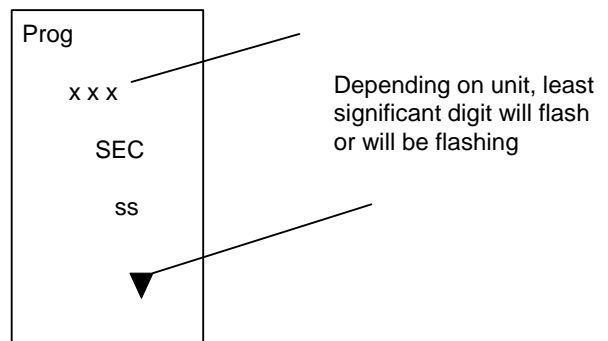


- E. To set delay time (next significant digit)

- Press "Prog" button until the desired # is set, when set press "Reset" button and the display will move to the next significant digit

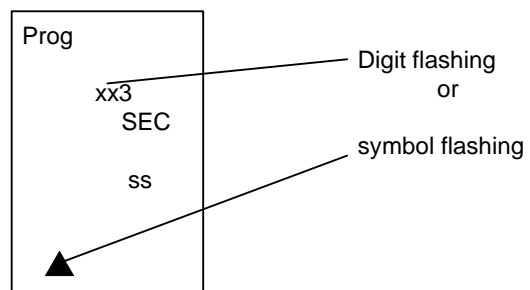


- F. To set delay time (most significant digit)
- Press "Prog" button until the desired # is set
 - Press "Reset" button

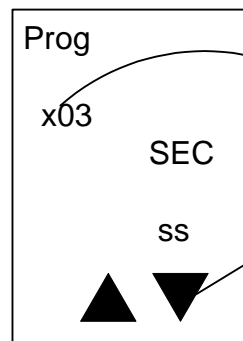


Digit flashing
or
Symbols flashing

- G. Set trigger type {(Level (LEU) or Edge (Edg))}
- Press "Prog" button to display Edg
 - Press "Reset" button when desired trigger type is indicated

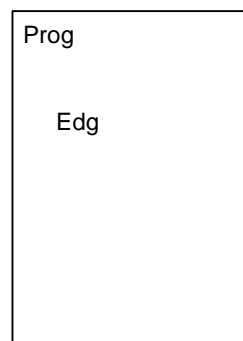


H. The timer now reviews your selected setting (if changes were made)

If the selection is correct
- Press "Reset" button

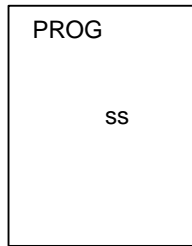


{Timer relay conf.} If incorrect press "Prog" button until you have the correct selection then press "Reset" {time scale}

If selection is correct
- Press "Reset" button

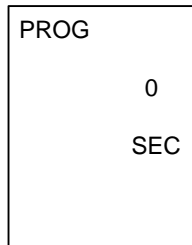


If incorrect press "Prog" button until you have correct selection then press "Reset" {time setting (lsd)}

If selection is correct
- Press "Reset" button

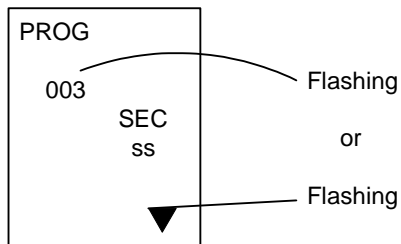


If selection is correct
- Press "Reset" button

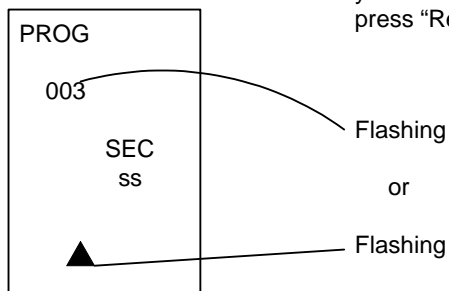

If incorrect press "Prog" button until you have correct selection - then press "Reset"


If selection is correct
- Press "Reset" button

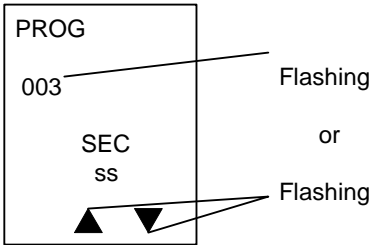

[Time Setting MSD]


If incorrect, press "Prog" button until
you have the correct selection, then
press "Reset"

If selection is correct
- Press "Reset" button

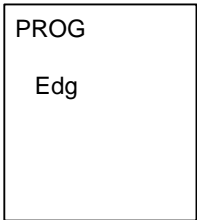

[Trigger type]

If incorrect, press "Prog" button until
you have the correct selection, then
press "Reset"

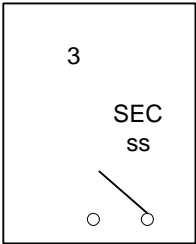

Operational state

(Rack active or not)